


Imaging Summary: TL_20260220_110k

Date	20260220
Folder	TL_20260220_110k
Start Time	20260220 12:20:51
End Time	20260220 14:04:33
Total Time (hrs)	1.73
Grid Squares Collected	2
Total Movies	332
Average Movies per Grid Square	166.0
Movies per Hour	191.91
Microscope	DFCI Tundra
Acceleration Voltage (kV)	100
Extractor Voltage (V)	3850
Spherical Aberration (mm)	1.6
Gun Lens	4
Spot Size	4
Intensity	0.448
EPU Version	3.15.0.11609
C2 Aperture (μm)	50
Objective Aperture (μm)	None
Camera	Ceta-F
Image Dimensions (pixels)	4096 x 4096
Nominal Magnification	110000
EPU Pixel Size (A/pix)	1.209
Calibrated Pixel Size (A/pix)	1.2
Beam Size (μm)	1.3
Pixel and Beam Size Calibration Date	20260110
Exposure Time (s)	2.864
Approx. Total Dose (e/pix)	43.09
Approx. Total Dose (e/A ²)	29.46
Approx. Dose Rate (e/pix/s)	15.05
Grid Type	HoleyCarbon
Grid Geometry	Square
EPU Measured Hole Size (μm)	1.32
EPU Measured Hole Center-to-Center Distance (μm)	2.6
Best Guess Hole Size and Spacing (μm)	1.2/1.3
Number of Acquisition Areas (Shots Per Hole)	1
AFIS	Yes
AFIS Clustering Distance (μm)	6.0
Number of Fractions	39
Defocus Values (μm)	[-2.1, -1.9, -1.7, -1.5, -1.3, -1.1, -0.9, -0.7, -0.5]


Notes

- Please contact Talya if any of these numbers appear to be incorrect! The script may need updating.
- The dose is approximated from the first movie. The total dose on specimen is slightly higher; however, if you did not record the dose when setting up collection, this will be appropriate for most (if not all) processing.
- The hole size and spacing is guessed based on the measure hole size function in EPU. If you are using an uncommon hole size/spacing, it may misidentify it.
- Pixel size is listed both as the pixel size automatically coded in EPU as well as the experimentally-calibrated pixel size. I advise that you use the calibrated pixel size in processing.



Atlas

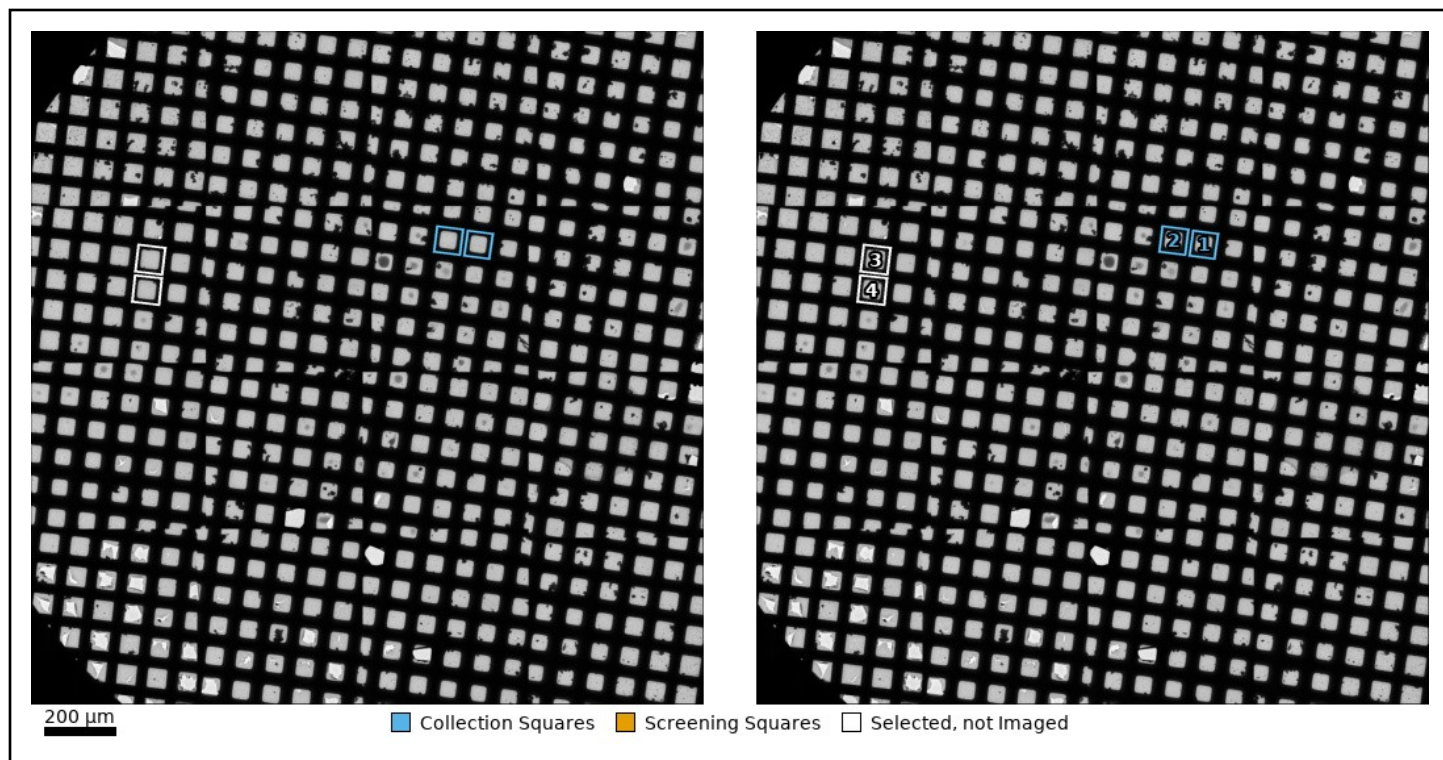

Template Definition

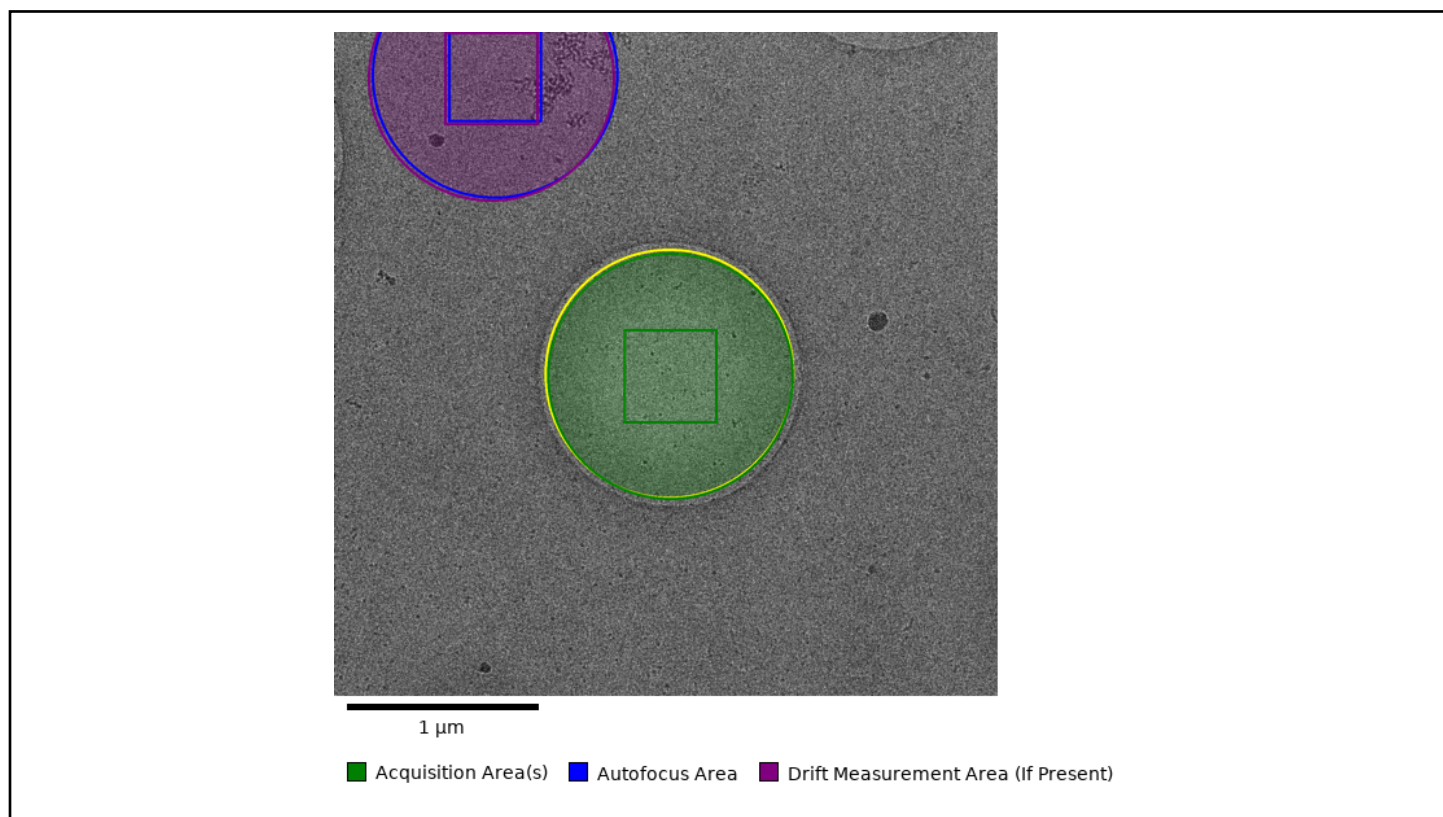


Grid Square 1 (EPU 19870688)

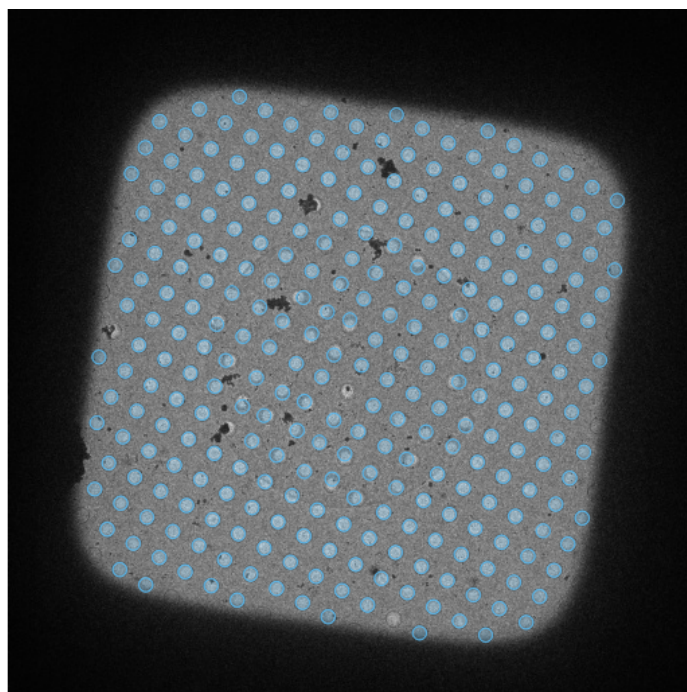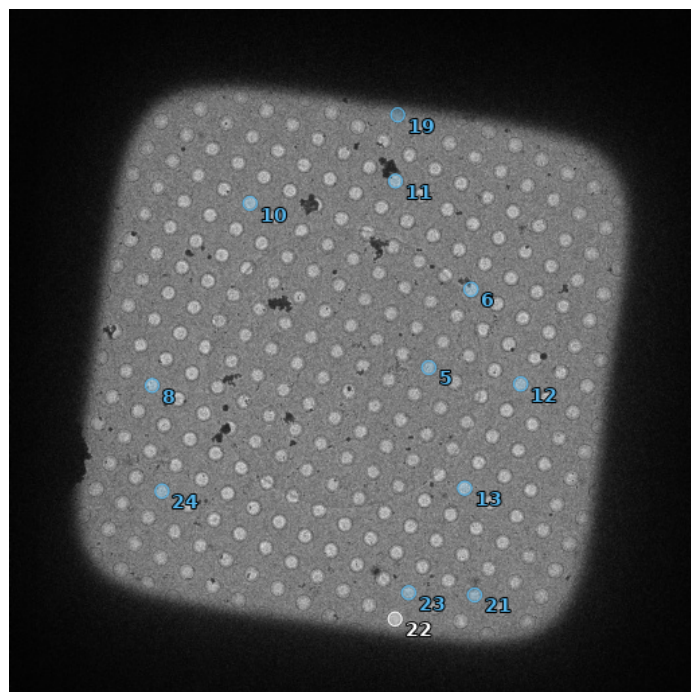


■ Collection Hole
 ■ Screening Hole
 ■ Selected, No Micrographs

A maximum of 12 holes (priority: screening, collection, selected no micrographs) are chosen for numbering and display

Foil Hole 1 (EPU 20017058)

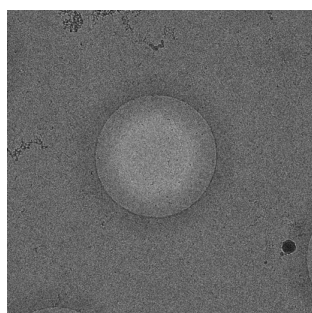


1 μm

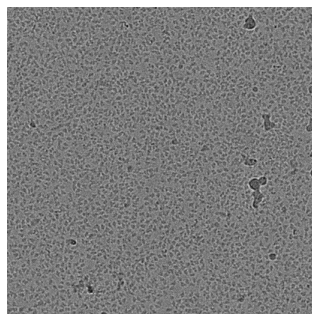


50 nm

Defocus: -2.10 μm

Foil Hole 2 (EPU 20017070)

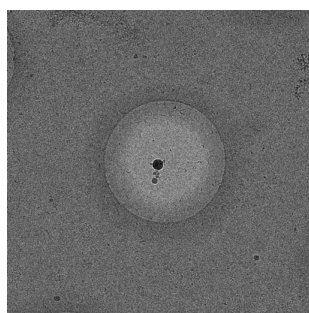


1 μm

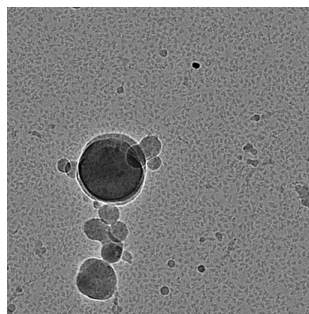


50 nm

Defocus: -1.90 μm

Foil Hole 3 (EPU 20017074)

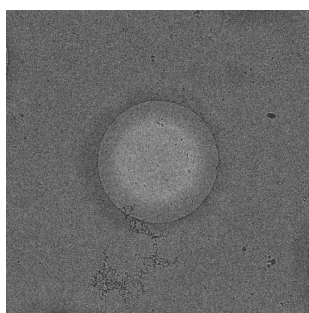


1 μm

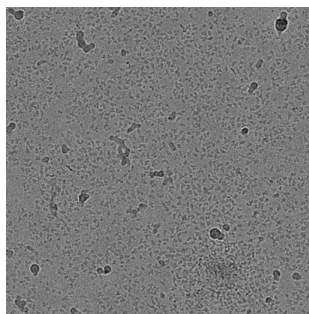


50 nm

Defocus: -2.10 μm

Foil Hole 4 (EPU 20017078)

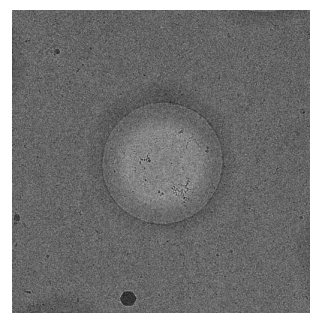


1 μm

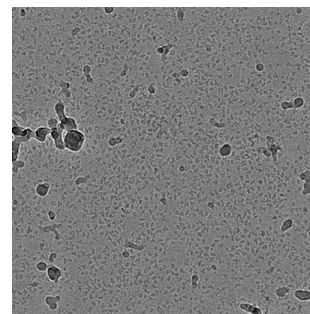


50 nm

Defocus: -1.90 μm

Grid Square 1 (EPU 19870688) (cont.)

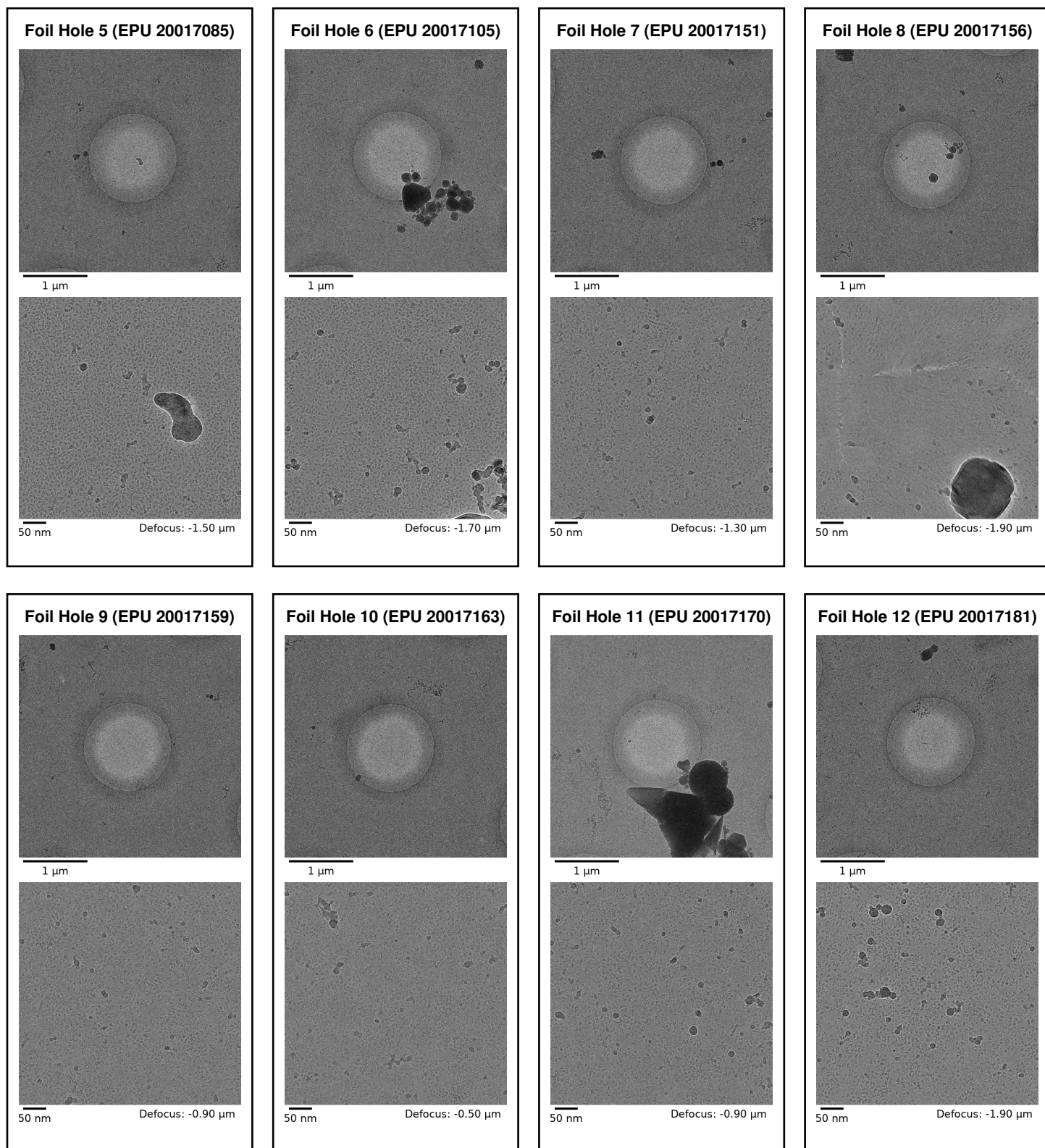


Grid Square 1 (EPU 19870688) (cont.)

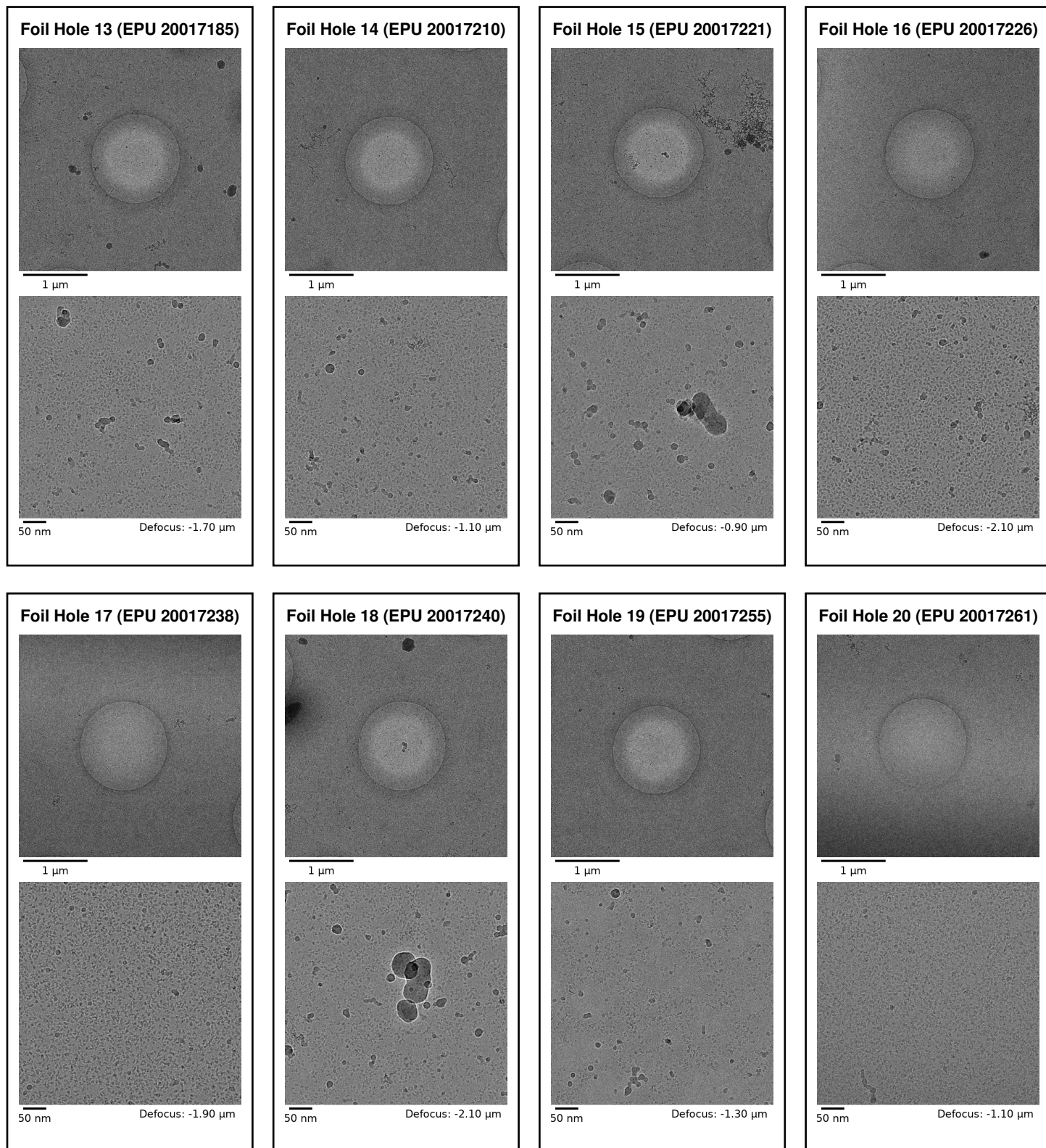


Grid Square 1 (EPU 19870688) (cont.)

Foil Hole 21 (EPU 20017273)

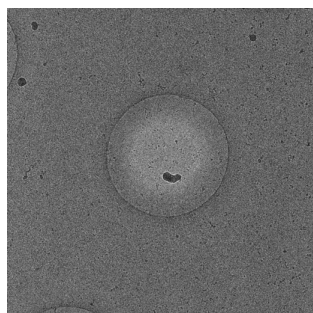


1 μm

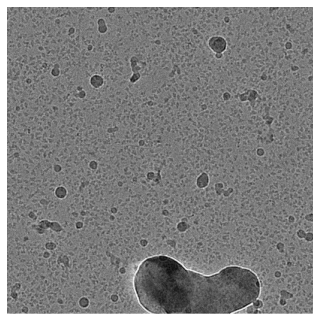


Defocus: -1.70 μm

Foil Hole 22 (EPU 20017279)

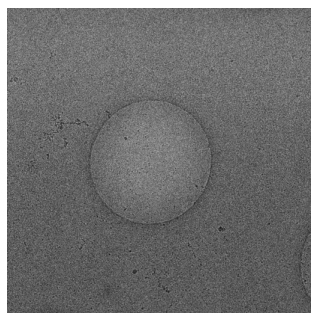

1 μm

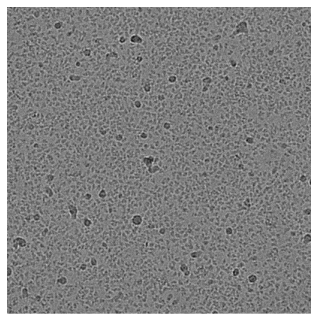


Defocus: -2.10 μm

Foil Hole 23 (EPU 20017282)




1 μm

Foil Hole 24 (EPU 20017291)

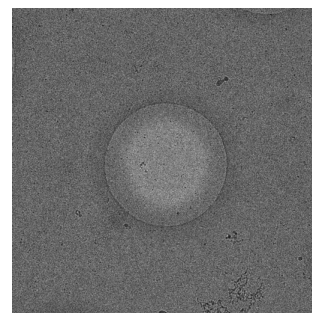


1 μm

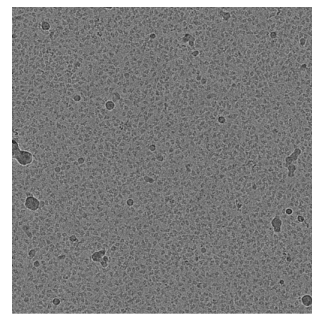


Defocus: -1.10 μm

Foil Hole 25 (EPU 20017299)

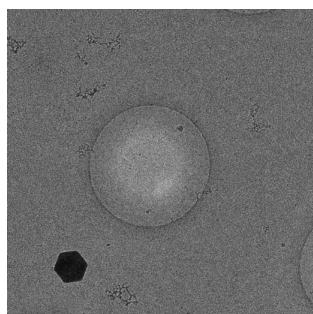


1 μm

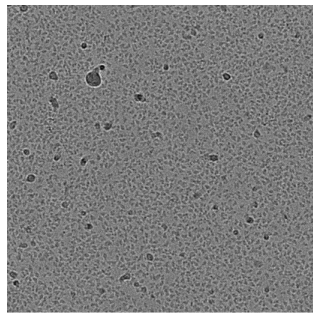


Defocus: -2.10 μm

Foil Hole 26 (EPU 20017303)

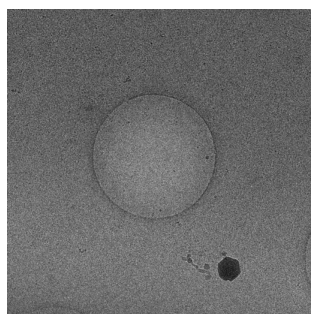


1 μm

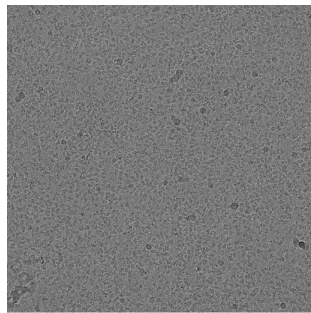


Defocus: -0.50 μm

Foil Hole 27 (EPU 20017307)

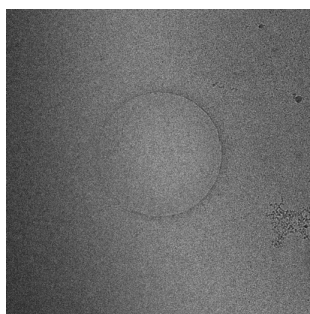


1 μm

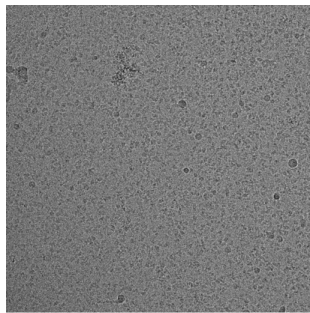


Defocus: -1.10 μm

Foil Hole 28 (EPU 20017312)

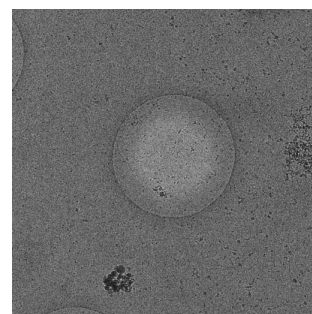


1 μm

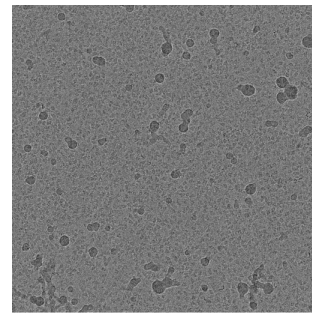


Defocus: -0.50 μm

Grid Square 1 (EPU 19870688) (cont.)

Foil Hole 29 (EPU 20017321)

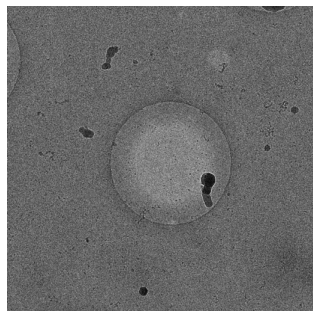


1 μm

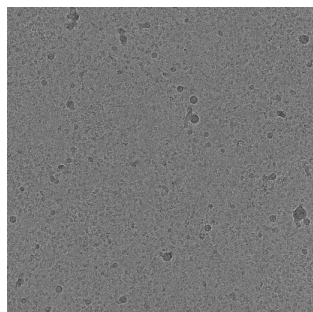


50 nm

Defocus: -0.50 μm

Foil Hole 30 (EPU 20017335)

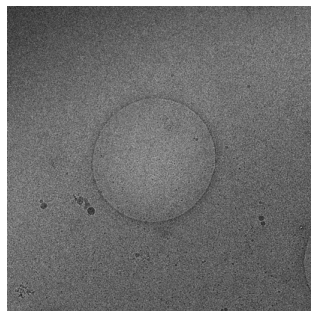


1 μm

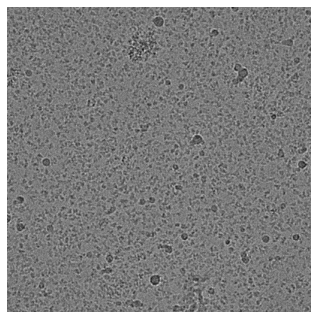


50 nm

Defocus: -1.90 μm

Foil Hole 31 (EPU 20017340)

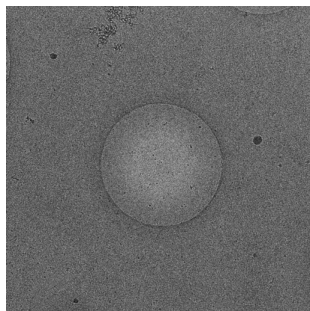



1 μm

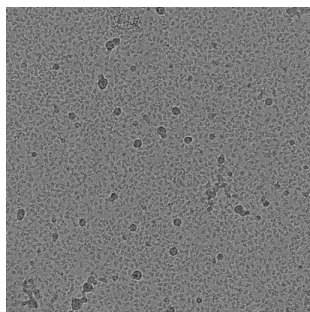

50 nm

Defocus: -1.30 μm


Grid Square 2 (EPU 19870686)

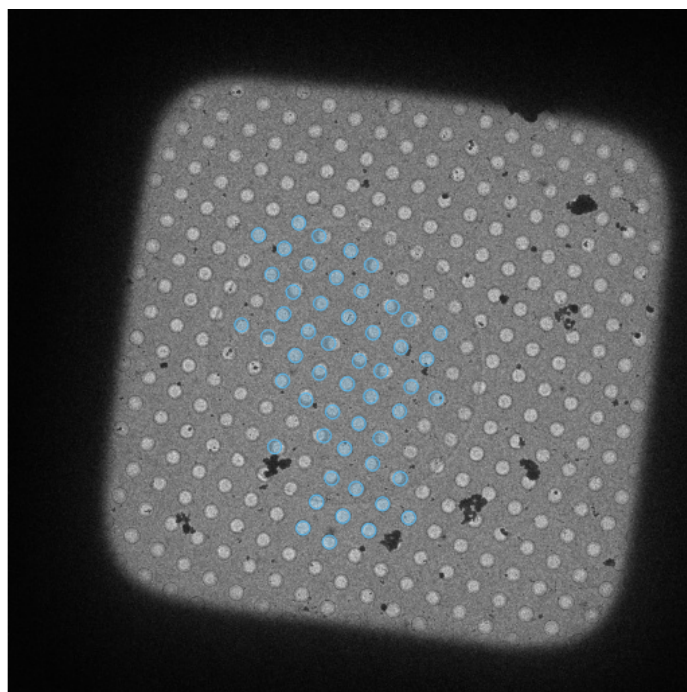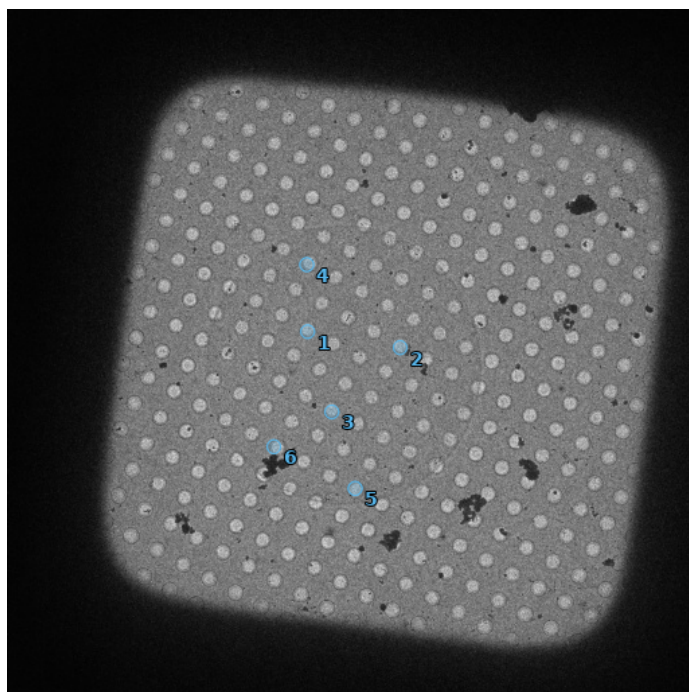


■ Collection Hole
 ■ Screening Hole
 Selected, No Micrographs

A maximum of 12 holes (priority: screening, collection, selected no micrographs) are chosen for numbering and display

Foil Hole 1 (EPU 20017661)

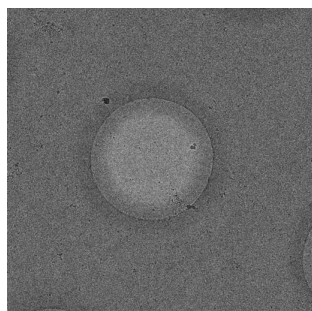


1 μm

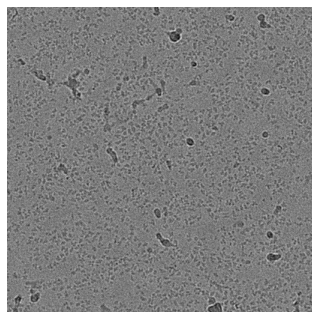


50 nm

Defocus: -2.10 μm

Foil Hole 2 (EPU 20017665)

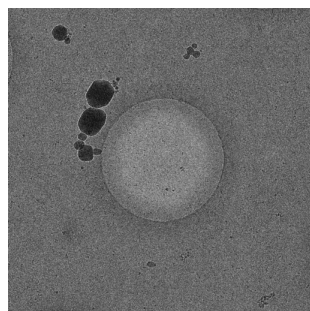


1 μm

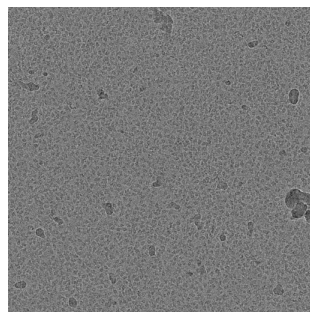


50 nm

Defocus: -0.50 μm

Foil Hole 3 (EPU 20017669)

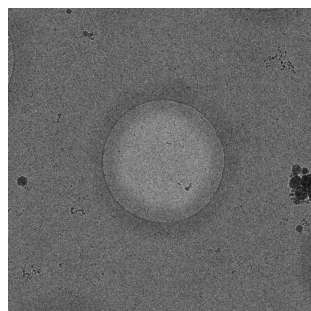


1 μm

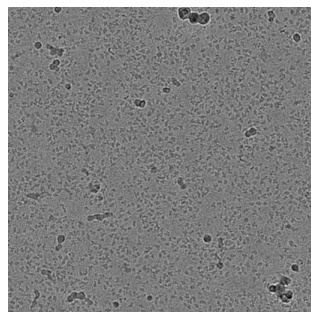


50 nm

Defocus: -2.10 μm

Foil Hole 4 (EPU 20017676)

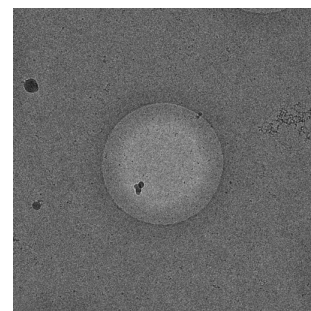

1 μm

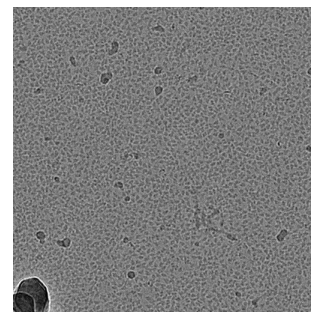


50 nm


Defocus: -1.90 μm

Grid Square 2 (EPU 19870686) (cont.)

Foil Hole 5 (EPU 20017713)

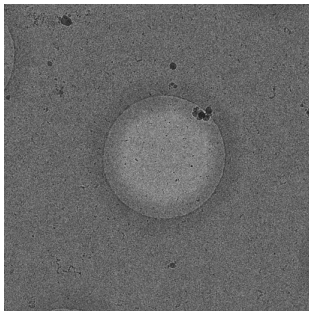


1 μm

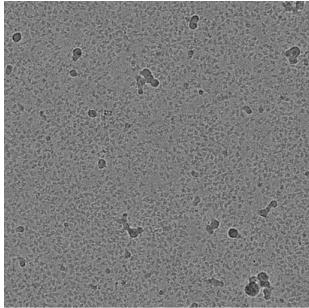

50 nm


Defocus: -1.50 μm

Foil Hole 6 (EPU 20017716)

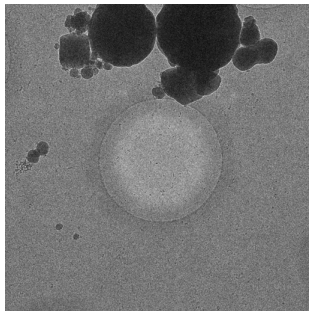

1 μm


Grid Square 3 (EPU 19870477)

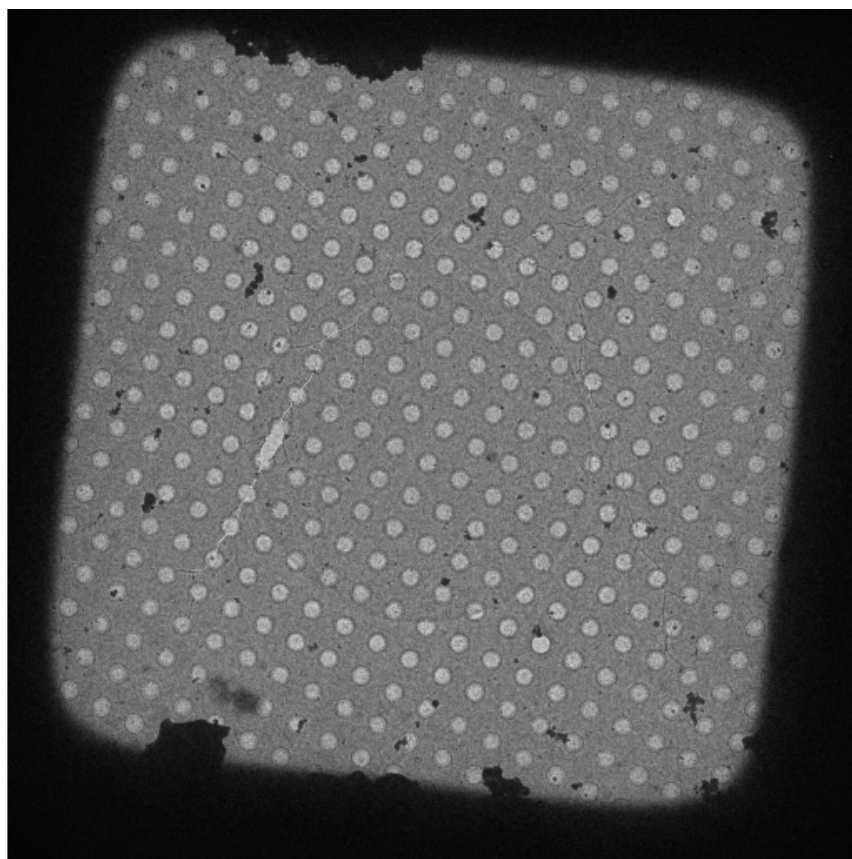

☐ Collection Hole ☐ Screening Hole ☐ Selected, No Micrographs

A maximum of 12 holes (priority: screening, selected no micrographs) are chosen for numbering and display

Grid Square 4 (EPU 19870491)

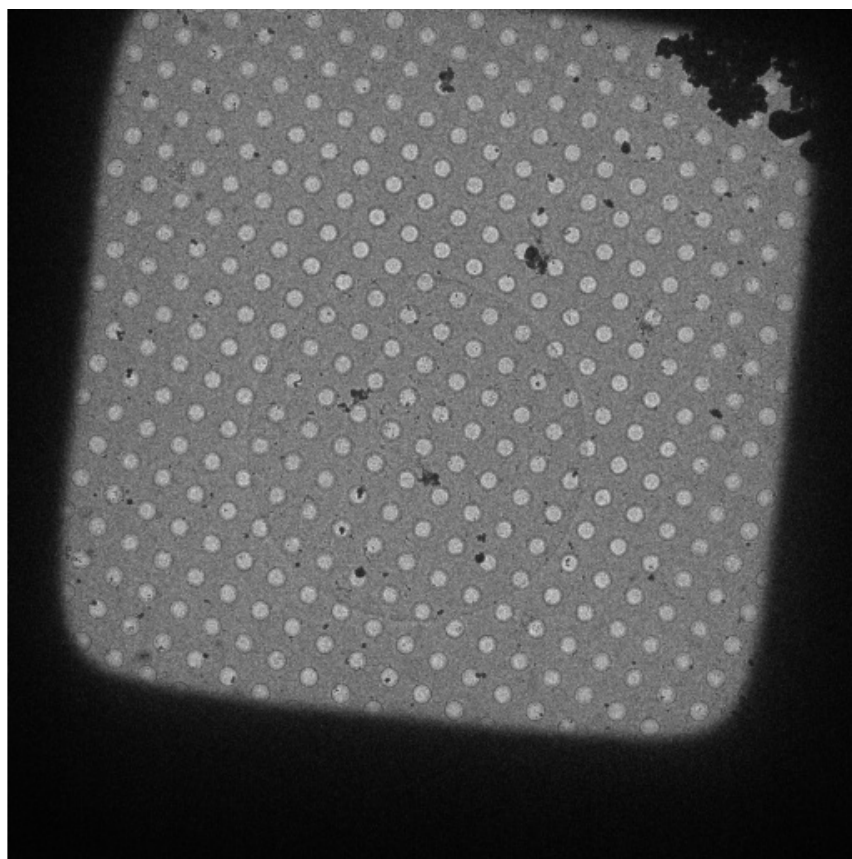

■ Collection Hole ■ Screening Hole □ Selected, No Micrographs

A maximum of 12 holes (priority: screening, selected no micrographs) are chosen for numbering and display